Dr. Mayer

LIFE & HEALTH

HOME DENTAL CENTER

Instruction Manual

Model: HDC5000

ENG
RO
BG
HU

IMPORTANT SAFEGUARDS

Please read the safety instructions before using this product. The safety instructions provided here help you use this product safely and correctly, and also to prevent you from injuring yourself or others.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not place or store the product where it can fall or be pulled into bathtub, sink or toilet.
- Do not place or drop the charger into water or other liquid.
- Do not reach for the charger base that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.
- Unplug the charger base before cleaning and ensure that it is dry before plugging it into an electrical outlet.

DANGER-VOLTAGE

- Do not plug the charger base into an outlet with a voltage other than specified on the bottom of the charger (100-240V).
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.
- This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

MEDICAL-WARNINGS

- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using the product.
- Consult your dentist if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).

WARNINGS


To reduce the risk of burns, electrocution, fire, or injury to persons:

- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the device by an adult responsible for their safety.
- Children should be supervised to ensure they do not play with the product.
- Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), please thoroughly clean your brush head and handle to avoid cracking the plastic.
- Certain liquid mouthwash may damage the device by causing the housing or tank to crack, when using mouthwash, rinse the tank and main unit with water immediately after use. Be especially careful when washing the hose connections at the main unit and the tank.
- Never fill the tank with hot water! The gum lines can be very sensitive and can be easily damaged by hot water.
- Use the unit within an ambient temperature range of 32F-104F.
- Never operate this product if it has a damaged cord, plug or charger. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the product away from heated surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact your dentist.
- Do not clean the product in the dishwasher.
- Consult your physician prior to using the toothbrush if you have any medical concerns.
- The toothbrush has been tested and is certified as per ETL regulations.
- Never drop or insert any object into any opening.
- For household use only. Do not use outdoors.

SAVE THESE INSTRUCTIONS

Water Flosser

MAIN UNIT

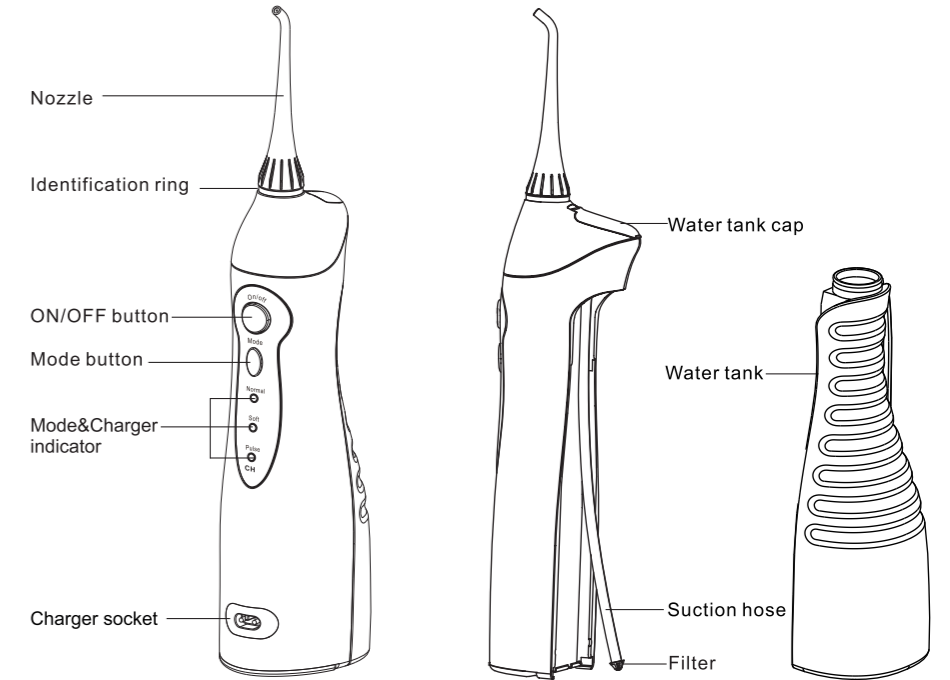ACCESSORY

GETTING STARTED

Charging the device (Fig. 1)

1. Insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug). (As ① on Fig. 1)
2. Insert the charging adaptor securely into the power outlet. (As ② on Fig. 1)
3. The "CH" indicator blinks when the device is charging, and it will be steadily on after the device is full charged.

Note:

If the "CH" indicator is blinking during or after an irrigating, that means the battery need to be recharged. An full recharging will take about 8 hours.

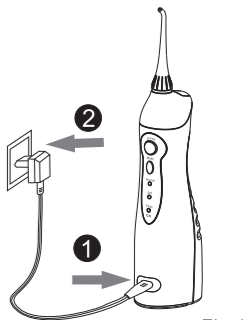
Fig.1

HOW TO USE

1. Attach nozzle onto the unit. (Fig. 2)
2. Open the water tank cap, hold the unit upright, fill the tank with water which temperature not exceed 104F, then close the cap.
3. Or, take down the water tank from the unit, hold the tank upright, fill the tank with water which temperature not exceed 104F, then attach the tank onto the unit. (Fig. 3)
4. Press "Mode" button to select the mode you desired. During irrigating, you can change the working mode.

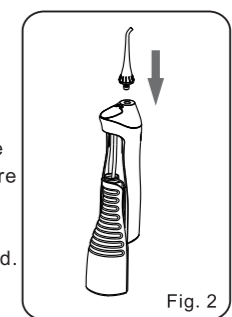
Fig. 2

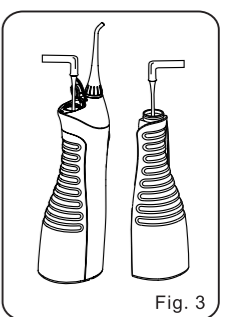
Fig. 3



MODE AND INDICATOR (Fig. 4)

Mode: Pressing the mode button will change the water jet setting and the respective mode indicator will light up.

Normal: for normal cleaning.

Soft: for gentle cleaning of sensitive gums.

Pulse: for cleaning and massaging the gums.

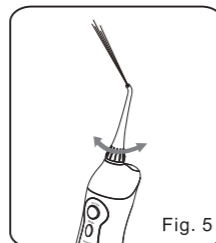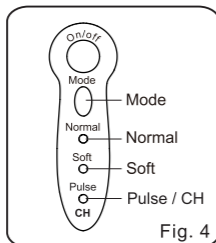5. Hold the unit in an upright position with nozzle tip point toward your teeth.

6. Keep your mouth opened slightly so that the water can flow out.

7. Press "ON/OFF" button to start irrigating. (Fig. 5)

- Position the water jet flow at a 90-degree angle to the teeth and gums.

- Move the water jet slowly along the teeth.

- Position the water jet along the gum line with the nozzle tip at a 90-degree angle to the gum tissue.

- Now clean between the teeth and around any braces, crowns or bridge-work thoroughly.

8. Press "ON/OFF" button again to turn off the device once irrigating complete.

9. Use of the Orthodontic Tip (Fig.6)

- Start on the last molar (back tooth). Aim water stream at 90° to the tooth at the gumline. Gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth. Continue process until you have cleaned the inside and outside of both the upper and lower teeth.

- Do not scrub or press hard, simply let the brush do all the work.

10. Use of the Tongue Cleaner (Fig.7)

- To use the tongue cleaner, set the unit to the lowest pressure setting. Gently place the tongue cleaner in the center/midline of your tongue, about halfway back and start the water flosser. With light pressure pull forward. You can increase pressure as you prefer.

AFTER USE

1. Open the water tank cap and dispose of the water in the tank. (Fig.8)





2. Press "On/Off" key to turn on the unit, and drain the water inside the appliance. (Fig.9)

3. Press "On/Off" key again to turn off the unit, and close the water tank cap.

4. Wipe the device with dry cloth.

Note:

To prevent bacteria propagating, make sure there is no remaining water inside the unit. If you're not going to use this appliance for long time, be sure wipe and dry the appliance ASAP before store it.

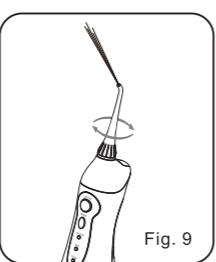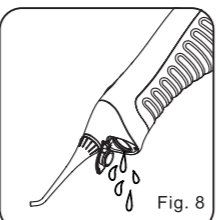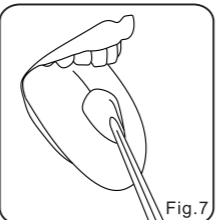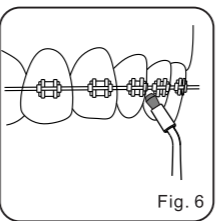MAINTENANCE

WATER TANK

Cleaning the device

Switch off the device and remove the water tank and the jet nozzle.

- Do not wash with hot water exceeding 122F.

- Use only water or a neutral cleaning agent to clean the main unit. Never use caustic or abrasive agents (e.g., vinegar or limescale removers), since they may damage the device.

MAIN UNIT

Use a cloth to clean the main unit. Wash hard-to-remove stains with water.

- Do not submerge the main unit in water.

- After washing with water, wipe dry with a cloth.

WATER TANK

Wash the water tank in water.

- Pour out the water inside the tank after cleaning.

- If you do not intend to use the device for a week or more, do dry the inside and outside of the tank.

- If cleaning requires a detergent, use a mild, neutral cleanser.

- Do not use harsh chemicals such as thinner, benzene or alcohol as these will discolor the surface. About once every six months wipe the prongs of the plug with a dry cloth.

TROUBLESHOOTING

Problem	Possible Cause	Solution
It does not work	No batteries inserted	Insert batteries in the unit
	The batteries in wrong direction	Re-insert batteries in correct direction
	The batteries are exhausted	Replace with new batteries
Water flow is very low	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	
	There is dirt stuck in the filter	Clean the filter
Water does not come out	The tank is empty	Fill up tank
	The appliance is tilted too much while in use	Use with the appliance in a vertical state

SONIC TOOTHBRUSH

BRUSH HANDLE

Power / Mode button

Mode indicator

Charging indicator

ACCESSORY

Professional Brush Head

Interdental Brush Head

CHARGER

Storage cover

Storage

Brush head stand

Charging cradle

Power cord

Power plug

GETTING STARTED

Charging the toothbrush

- Be sure to turn off the toothbrush before charging. (Fig.1) Place the handle in the charger. The charging indicator will light up blue. (Fig.2)

Notes:

- It takes at least 24 hours to charge the battery fully.



- The toothbrush must be charged for 24 hours, if:

- this is the first time using the toothbrush;
- the toothbrush has not been used for more than 3 months;

- To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use. It cannot overcharge.

- During charging, the toothbrush will get warm, but this is normal.

- The operating duration will become shorter as the battery gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.

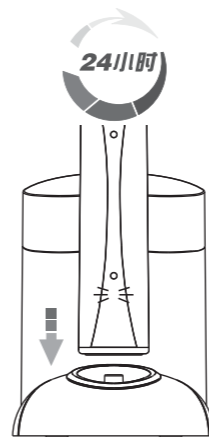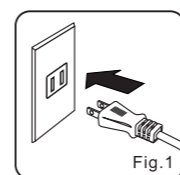HOW TO USE

1 Slide the brush head onto the handle. Push brush head down until it locks into place on handle.

Notes: There is small gap between color ring and the handle. (Fig. 3)

2 Wet the bristles and apply a small amount of toothpaste.

3 Place the bristles against the teeth at a slight angle towards the gum line. (Fig. 4)

4 Press the Power / Mode button to turn on the toothbrush, It starts working at Normal mode.

Normal (H): for normal cleaning.

Soft (S): for gentle cleaning of sensitive gums.

Massage (M): for cleaning and massaging the gums.

Notes: During brushing, if press the Power / Mode button, toothbrush will work at next mode and start another 2 minutes cycle. The function turn of the button is: power on at normal mode, soft mode, massage mode and power off.

5 Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.

6 In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use. (Fig. 5)

7 Use of the Interdental Brush Head

- Specially designed shape is ideal for cleaning braces, crowns or spaces between poorly aligned teeth.

- Place the brush lightly and vertically against your teeth. Guide the brush head move slowly along the gap between the teeth. (Fig. 6)

8 Pull brush head to detach from the toothbrush.

9 Rinse both the brush head and handle under running water wipe them dry. (Fig. 7)

10 Put the handle back into the charger.

MAINTENANCE

Be sure to remove the plug from the electrical outlet before cleaning the product.

Charger

- Pull the storage to remove. (Fig. 8)

- Using a soft cloth, wipe off any dirt on the charger and the storage. (Fig. 9)

- Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

Toothbrush

- Remove the brush head.

- Wipe all surface areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.

- Rinse the brush head and bristles after each use.

- Do not clean in the dishwasher.

- Any other servicing should be performed by an authorized service representative or that the product has no user serviceable parts.

TROUBLESHOOTING


Problem	Possible cause	Solution
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 24 hours.

DISPOSAL

The power source of this product is nickel-metal hydride battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.

Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions. The appliance must be disconnected from the supply mains when removing the battery.

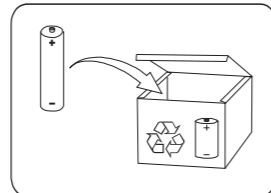Water Flosser

1. Remove 3pcs rubber feet on the bottom the handle.
2. Insert a screwdriver into the slot located at the bottom of the handle and turn counterclockwise to release the bottom cap.
3. Pull out the batteries from the handle.

Sonic Toothbrush

- Before removing the battery, be sure to completely drain the unit of power.
- Use a screwdriver to remove the screw and battery cover. (Fig.10)
- Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle. (Fig.11)
- Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier. (Fig.12)
- Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.
- Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

Fig.10

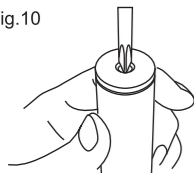
Fig.11



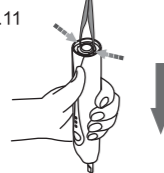
Fig.12

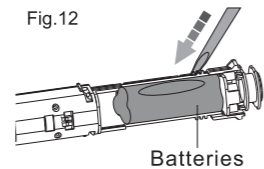DENTAL MIRROR

Use of the Dental Mirror

- 1 Allows vision in certain locations of the mouth where visibility is difficult or impossible.
- 2 By pressing the power button the LED mirror light will light up. Press again, the LED mirror light will go off.

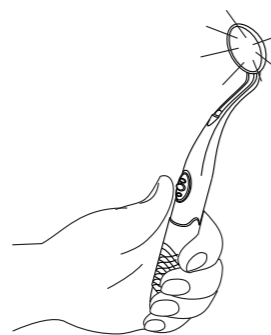Dr Mayer
Life and Health
E-mail: office@doctor-mayer.com
Web: www.doctor-mayer.com

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